97050590
rev. 001
10/2015

ITALIANO (Lingua originale).....	3
ENGLISH (translation from the original in Italian).....	17
FRANÇAIS (traduction de l'original en italien).....	31
DEUTSCH (Übersetzung aus dem Original in Italienisch).....	45
ESPAÑOL (traducción del original en italiano).....	59
PORTUGUÊS (tradução do original em italiano).....	73
ΕΛΛΗΝΙΚΑ (αση από το πρωτότυπο στα ιταλικά).....	87
РУССКИЙ (Перевод с оригинала в итальянском).....	101
POLSKI (tłumaczenie z oryginału w języku włoskim).....	115
TÜRKÇE (İtalyan orijinalinden çeviri).....	129

REF 97050590

C7/C8/C9

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INDICE

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CONTENTS

1. General information	18
1.1. Symbols	18
1.2. Shipment and storage conditions	18
1.3. Disposing of the packing materials	19
1.4. Disposing of the dentist's stool	19
1.5. Warranty	19
2. Main features	21
2.1 Use of the dentist's stool	22
2.1.1. Reference standards	22
2.1.2. Environmental conditions	22
2.2. Identification	22
2.3 Technical information	22
3 Assembling the dentist's stool	23
4. Using the dentist's stool	23
4.1. Dentist's stool, model C7	23
4.2. Dentist's stool, model C8	24
4.3. Dentist's stool, model C9	25
5. Maintenance	27
5.1. Cleaning and disinfecting the stool	27
5.2. Replacing the castor wheels	28
5.3. Replacing the seat upholstery	28
6. Trouble-shooting	29


1. General information

- This manual describes the C7-C8-C9 model dentist's stool. It contains important instructions regarding their operation, functions and maintenance as well as a trouble-shooting guide to any problems which might arise
- No part of this manual is to be reproduced, stored in a retrieval system or transmitted in any form or by any means, i.e. electronic, mechanical, photocopying, translation or otherwise, without the prior written permission of CEFLA s.c.
- All information, specifications and illustrations contained herein are subject to change.
As CEFLA s.c. has a company policy of continual development, it reserves the right to make changes without prior notice.
- The original text of this manual is in Italian.

1.1. Symbols

The following symbols and messages have been used in this manual:

- 1) WARNING/CAUTION! This indicates a situation where failure to observe the instructions may cause damage to the stool or injury to the user.
- 2) NOTE: This indicates important information for the user and/or service personnel.
- 3) Symbol indicating separate collection for electrical and electronic equipment in compliance with Directive 2012/19/UE.
- 4) Manufacturer.
- 5) Month and year of constructions.
- 6) Apparatus serial number.

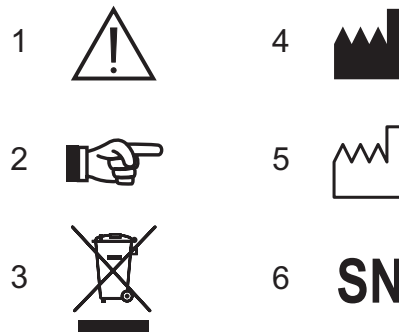1.2. Shipment and storage conditions

When packed for storage or shipment the dentist's stool can be kept under the following ambient conditions for no more than 15 weeks:

- Ambient temperature from - 10 to +70°C
- Relative humidity from 10 to 90% including condensation
- Atmospheric pressure from 500 to 1060 hPa (from 500 mbar to 1060 mbar).


The shipment, handling and unpacking instructions are printed on the cardboard box.

WARNING!

Carefully follow these instructions.

- 10) Observe the arrows when handling the packing containers.
- 11) Avoid banging the boxes.
- 12) Protect from damp.
- 13) Do not use hooks when handling.
- 14) Observe nameplate indicating ambient conditions required for storage.

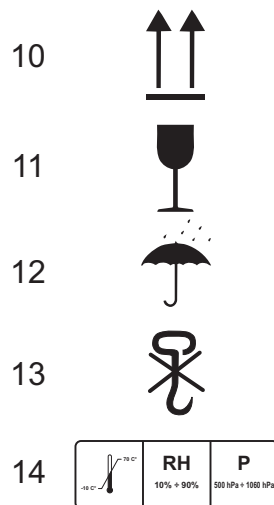1.3. Disposing of the packing materials

All the packing materials used are environ-mentally-friendly and 100% recyclable:

- cardboard,
- polythene bag,
- polyethylene jiffy paper
- polystyrene

Recycling the packing materials helps save raw materials and reduce waste.

Deliver the packing materials to an authorized recycling center.

1.4. Disposing of the dentist's stool

Disposing of the dentist's stool poses no threat to the environment.

If the stool is no longer to be used, it should be delivered to an authorized recycling center in order to protect the environment.

1.5. Warranty

CEFLA s.c. guarantees its products are constructed in a state-of-the-art manner, only top quality materials are used, all the required tests and inspections have been carried out and the product conforms to all current applicable standards.

CEFLA s.c. warrants products covered by this warranty to be free from defects in materials and workmanship.

WARRANTY GENERAL TERMS AND CONDITIONS

- 1) **COVERAGE:** the product is covered by a 12 (twelve) month warranty from the date of installation at the customer's location. In the event that the device works together with CEFLA s.c.-made dental equipment, refer to the dental equipment warranty conditions to know warranty duration.
- 2) **EXCLUSIONS:** this warranty does not cover:
 - a) faults caused by failure to perform routine maintenance or unintended, inappropriate or misuse of the product;
 - b) periodic checks and maintenance;
 - c) repair or replacement of fragile parts, parts subject to wear and tear or whose service life cannot be estimated unless there is proof that they were defective at the time of delivery;
 - d) all labor, travel and lodging expenses;
 - e) failures or break downs caused by unintended, inappropriate or misuse;
 - f) failures or break downs due to extraordinary chemical or electrical phenomena;
 - g) failures or break downs resulting from use of cleaning, disinfecting or sterilization chemicals or processes not expressly indicated in the operating and maintenance manual;
 - h) normal wear and tear of plastic parts .
- 3) **LIMITS:** defective parts will be repaired or replaced free of charge. However, this warranty does not cover replacement of the entire unit. As far as components applied or installed made by CEFLA s.c. or third parties that carry their own warranty are concerned, the conditions, terms, limits and exclusions given on the certificates are valid.

The date the warranty expires will not change even if one or more repairs are made during the warranty period. The purchaser may not make any damage claims for equipment downtime.
- 4) **WARRANTY CLAIMS:** if the purchaser makes a warranty claim regarding the quality or conditions of the equipment received, the purchaser may not, in any case, stop and/or delay payment.
- 5) **DISPUTES:** any disputes concerning application and interpretation of the Warranty Certificate will be judged by the Court of Bologna (Italy), regardless of where the purchase contract was signed.
- 6) **EXCEPTIONS:** particular exceptions granted in special cases as regards this warranty do not imply any rights given to the Purchaser and are to be considered valid only for the specific case.
- 7) **OTHER:** Regarding any aspects not expressly covered in this Warranty Certificate the Civil Code and applicable Italian laws shall rule.

LOSS OF WARRANTY

The warranty will be void if:

- a) the equipment shows signs of damage due to being dropped, exposure to flames, liquid spills, lightning, natural disasters, atmospheric phenomena or any other causes that cannot be traced to defects in workmanship;
- b) the equipment was not installed in accordance to CEFLA s.c. installation instructions or installed by personnel not duly authorized;
- c) the equipment was serviced, repaired, altered or modified in any way by the Purchaser or third parties that have not been duly authorized;
- d) when the claim for warranty repair was made the serial number on the product has been removed, canceled, counterfeited, etc...
- e) the Purchaser stops and/or delays payment of any sum for the purchase and/or maintenance of the equipment for any reason whatsoever;
- f) scheduled maintenance or other operations specified in the operating and maintenance manual are not carried out as indicated.


2. Main features

This chapter gives an overall idea of the main features of the C7-C8-C9 model dentist's stools

①

C7 dentist's stool


- Wheeled base (5 castor wheels with brake)
- Base with footrest ring (optional).
- Symmetrical seat which can be used on either side of the patient's dental stool.
- Tiltable backrest which can be adjusted as required by using the knob provided.
- The backrest can also be adjusted in height (it can be moved up and down approx. 60 mm).

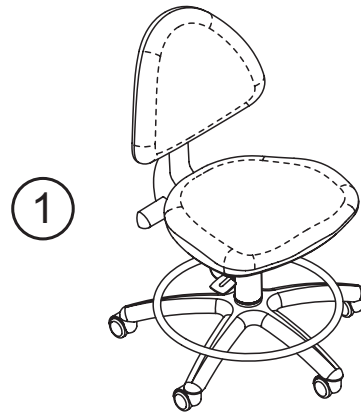②

C8 dentist's stool


- Wheeled base (5 castor wheels with brake)
- Base with footrest ring (optional).
- Symmetrical seat which can be used on either side of the patient's dental stool.
- Adjustable backrest (it can be moved to the left or right).

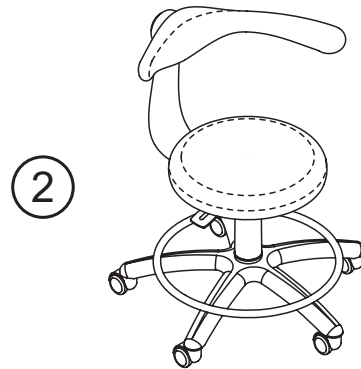③

C9 dentist's stool.

- Sliding base on 5 braked wheels.
- Base with footrest ring (optional).
- Perfectly symmetrical dentist's stool that can be used both to the right and left of the dental stool.
- Backrest horizontally adjustable to 5 positions according to build (total excursion 50mm);
- Backrest vertically adjustable to 5 positions (total excursion 60mm);
- Backrest 20° tiltable to ensure lumbar support;
- Seat front 20° tiltable to avoid pressure on the thighs.

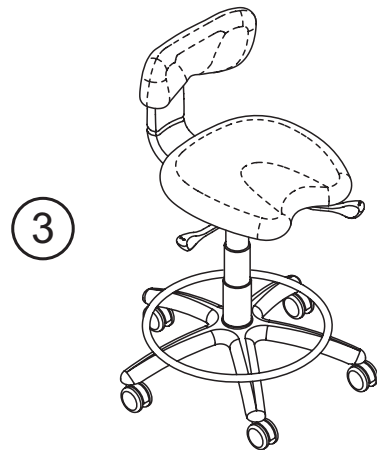2.1 Use of the dentist's stool

The C series dentist's stools are accessories for medical use designed to ensure that the dentist or the assistant is properly supported during normal dental practice

WARNING!

- The maximum allowable load that the stool can bear is 135 Kg. Do not exceed this limit.
- If the rear wheels are obstructed while moving back, the stool may overturn.
- The device is not intended to be used in rooms subject to the risk of explosion.

2.1.1. Reference standards

Reference standards: The C series dentist's stools are devices designed in compliance with the ISO 7493 standards.

2.1.2. Environmental conditions

The C series dentist's stools must be used in environments with the following characteristics:


- Temperature between 10 and 40°C;
- Relative humidity between 30 and 75%.

2.2. Identification

The identification plate is positioned underneath the seat..

Data plate information:

- Manufacturer.
- Stool model.
- Serial number.
- Month and year of construction.

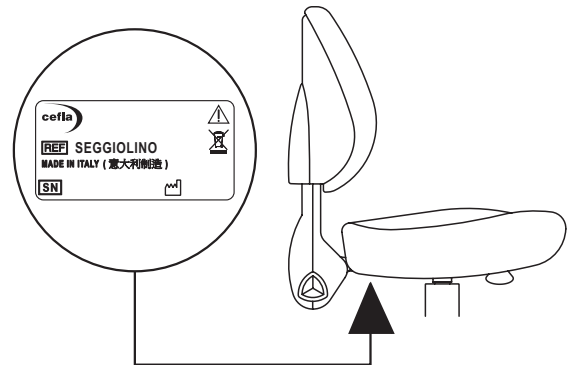2.3 Technical information

Dentist's stool, model:	C7	C8	C9
Base wheels outside diameter (mm):	605	605	605
Seat height from the floor (mm):	420	450	490
Seat vertical adjustment travel (mm):	130	140	165
Backrest horizontal adjustment travel (mm):	-	-	50
Backrest vertical adjustment travel (mm):	60	-	60
Seat dimensions (mm):	460x420	ø380	410x420
Backrest adjustment angle:	+11° / -9°	-	+10° / -5°
Seat front adjustment angle:	-	-	0° / -20°
Gross weight (kg):	18	18	18
Net weight (kg):	16	16	16
Packaging dimensions (mm):	570x570x990(h)	570x570x990(h)	570x570x990(h)

3 Assembling the dentist's stool

The C7-C8-C9 dentist's stools are supplied in two parts inside the packing container.


The following items should be found inside the packing container:

- N°1 dentist's stool.
- N°1 User's manual.

Steps to be taken to assemble the stool:

- Take the bottom part (A) of the stool out of the packing container and place it on the floor.
- Take the top part (B) of the stool out of the packing container and install it on the shank of the gas-operated spring.
- Secure the two parts together by pressing firmly down on the seat (e.g. actually sit on it).

When the above steps have been completed, the stool is ready for use.

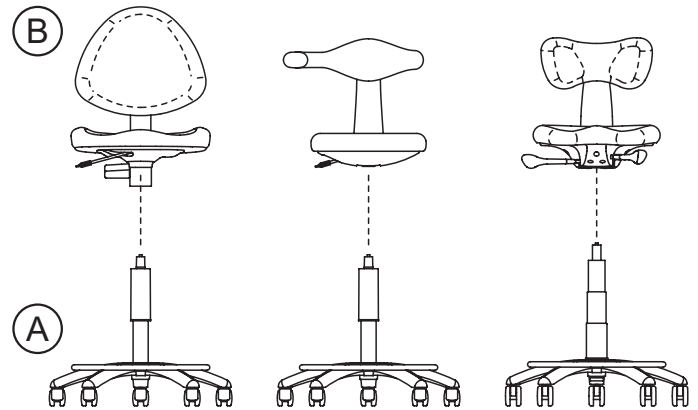4. Using the dentist's stool

Before reading this chapter it would be a good idea to identify exactly which model dentist's chair you have purchased - see Paragraph 2. To avoid confusion, only read the paragraphs relevant to the stool model you have purchased.

4.1. Dentist's stool, model C7


Adjustment levers.

S Knob used to adjust the backrest tilt.

L Seat up/down adjustment lever.


Seat height adjustment.

- Sit on the seat as the weight of your body will help move the seat down.
- Keep the lever (L) pulled up and the seat will move down with the help of your body weight.
- To move the seat up again, keep the lever (L) pulled up and lighten the load on the seat by standing up slightly.

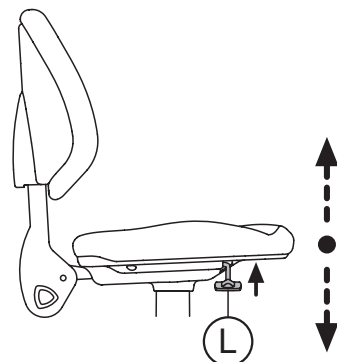Adjusting the backrest tilt.

- Lean away from the backrest to remove your weight from it.
- Turn knob S clockwise to tilt the backrest forwards.
- Turn knob S anti-clockwise to tilt the backrest backwards.

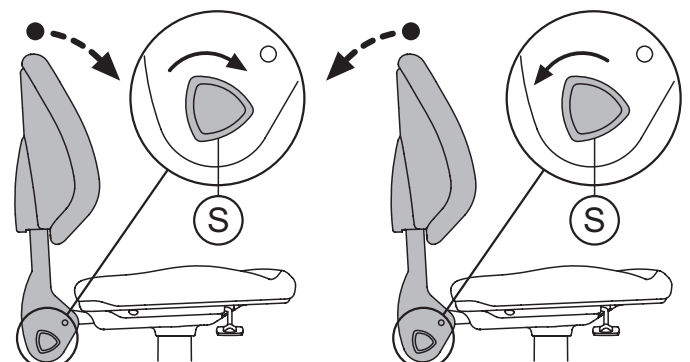Adjusting the height of the backrest

The height of the backrest can be adjusted (it can be moved a maximum of 60mm.).


The adjustment can be made manually by pulling the backrest up.

NOTE: There are nine positions the backrest can be set in.

To lower the backrest, proceed as follows:

- Move it up to the maximum height until the release mechanism is activated.
- Push the backrest down again until it is at the minimum height possible.
- Move the backrest to the required position.

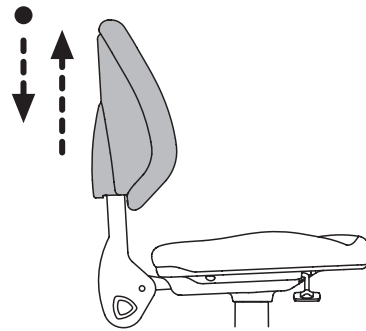
4.2. Dentist's stool, model C8

Adjustment levers.

L Lever to move the seat up and down.

Seat height adjustment.



- Sit on the seat as the weight of your body will help move the seat down.
- Keep the lever (**L**) pulled up and the seat will move down with the help of your body weight.
- To move the seat up again, keep the lever (**L**) pulled up and lighten the load on the seat by standing up slightly.

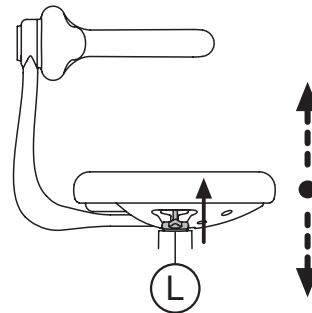Adjusting the backrest.

It is advisable to adjust the position of the backrest of the dentist's stool from behind.

- Simply move the backrest round on its support post in order to bring it into the required position.
- To change the backrest over from its RH position to the LH one, or viceversa, simply take hold of it with both hands near the support and push it forwards. In this way the lock mechanism disengages and the backrest can then be turned round to the new position where it will lock in place automatically.

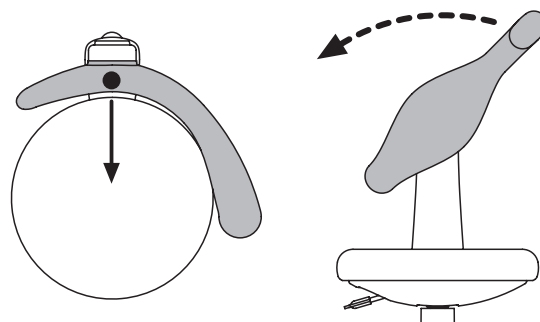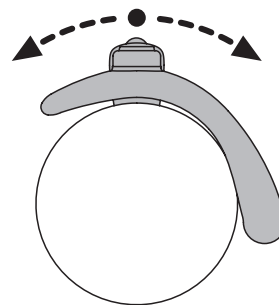
4.3. Dentist's stool, model C9

Adjustment levers.

- L** Seat up/down adjustment lever (can be operated with the back of the hand to avoid touching it with the fingers).
- X** Seat front tilt adjustment lever.
- S** Backrest horizontal adjustment lever.

Seat height adjustment.

- Sit down on the dentist's stool so that you can use your own body weight to help push it down.
- Push the lever (**L**) horizontally towards the rear and the seat will go down aided by your weight bearing down on it.
- To move the seat back up, hold the lever (**L**) pushed horizontally towards the rear, stand up and the seat will automatically move back up.

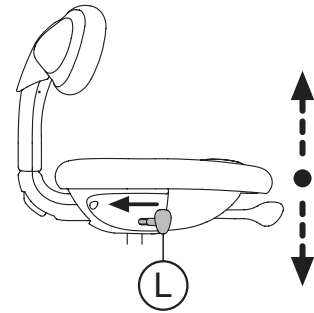
Seat front tilt adjustment.

This adjustment allows obtaining a perfect angle of 110° between the thigh and the leg:

- Push the lever (**X**) horizontally to adjust the seat front tilt.

 **NOTE:** Pushing the lever backward, the seat tilts downward (-20°); pushing the lever forward, the front part of the seat aligns with the rear part (0°).

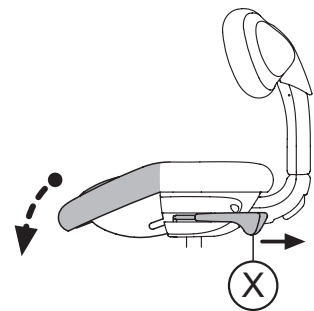
Backrest horizontal sliding adjustment.

This adjustment allows adapting the backrest position according to build.

- Pull the lever (**S**) upward and push or pull the backrest structure.


 **NOTE:** The backrest can be horizontally adjusted to 5 positions for a total travel of 50 mm.

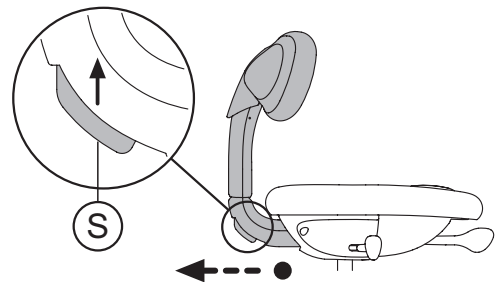Backrest tilt.

The backrest can be tilted on the vertical plane for lumbar support of the operator.

The backrest can be turned in a range between +10° and -5°.

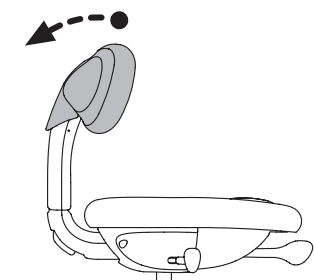Backrest vertical sliding adjustment.


This adjustment allows adapting the backrest position according to build:

- The adjustment is made by pulling the backrest upward.

To lower the backrest, operate as follows:

- Pull the backrest up to its maximum height to engage the release mechanism.
- Then push the backrest down until it reaches its minimum height.
- Position the backrest in the desired position.

 **NOTE:** The backrest can be vertically adjusted to 5 positions for a total travel of 60 mm.

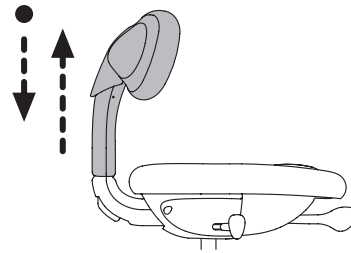Advice to properly adjust the dentist's stool.

Adjust the stool so that the lumbar support follows the lumbar lordosis of your spine (fig.3) without accentuating it (fig.2) nor slouching the back (fig.1).



To this end, the lumbar support should rest on the lumbar vertebrae close to the pelvis:

A: pelvis;

B: spine.

To adjust the stool to the most comfortable position, we suggest you:

- Sit down on the stool and first of all adjust the seat height using the right-hand lever (**L**);
- Then adjust the seat front tilt using the left-hand lever (**X**);
- Next adjust the lumbar support first in height and then, if necessary, horizontally with respect to the seat using the rear handle (**S**).

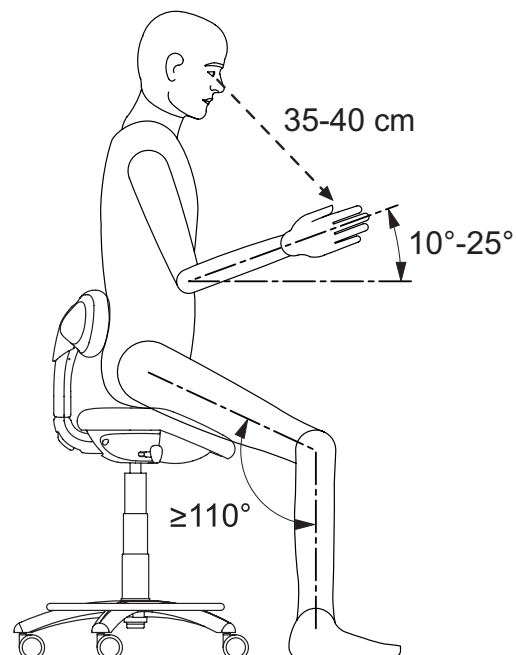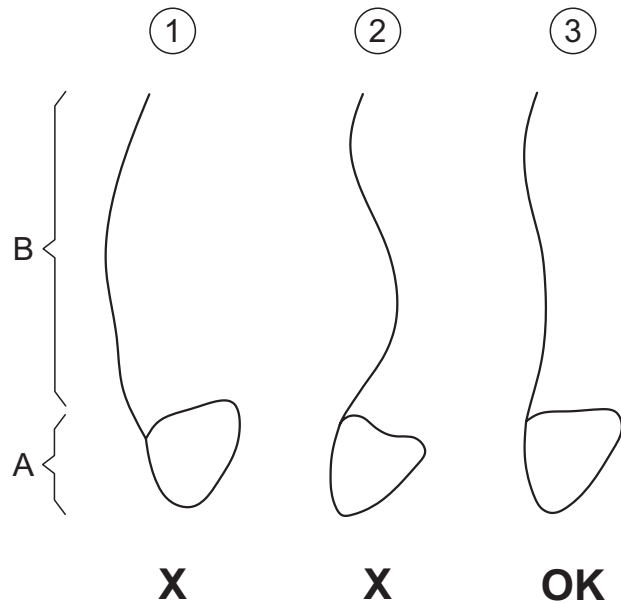5. Maintenance

The C7-C8-C9 dentist's stools have been designed so as not to require any maintenance during their useful life. No scheduled maintenance is therefore required. Nonetheless, it is advisable to have the dentist's stool checked by qualified CEFLA s.c. staff at least every two years.

The established useful life of the device is 10 years. After this period of time, the stool needs to undergo extraordinary maintenance in order to check and replace any worn components.

WARNING!

If the dentist's stool needs to be repaired, it is advisable to have this done by specialised personnel from our After-Sales Service Dept.

5.1. Cleaning and disinfecting the stool

Cleaning is the first step of any disinfecting process. Physically scrubbing with detergents and surface-active substances and rinsing with water removes a considerable amount of micro organisms. If a surface is not clean first, the disinfecting process cannot be successful. If a surface cannot be adequately cleaned, it should be protected with barriers.

The outer parts of the equipment must be cleaned and disinfected using a product for hospital use with indications for HIV, HBV and tuberculocide (medium-level disinfectant) specifically for small surfaces.

The various drugs and chemical products used in dentist's surgeries may damage the painted surfaces and the plastic parts.

Research and tests run show that the surfaces cannot be fully protected against the harsh action of all products available on the market. We therefore recommend protecting with barriers whenever possible. The harsh actions of chemical products also depend on the amount of time they are left on the surfaces. It is therefore important not to leave the product on the surfaces longer than the time specified by the manufacturer.

It is recommended to use the specific medium-level disinfectant, **STER 1 PLUS** (CEFLA s.c.), which is compatible with:

- **Coated surfaces and plastic parts.**
- **Upholstery.**
- **Uncoated metal surfaces.**

If you do not use STER 1 PLUS, it is recommended to use products that contain at maximum:

- **Ethanol 96%. Concentration: maximum 30 g per 100 g of disinfectant.**
- **1-propanol (N-propanol, propyl alcohol, N-propyl alcohol). Concentration: maximum 20 g per 100 g of disinfectant.**
- **Combination of ethanol and propanol. Concentration: the combination of the two should be maximum 40 g per 100 g of disinfectant.**

ATTENTION!

- **Do not use products containing isopropyl alcohol (2-propanol, iso-propanol).**
- **Do not use products that contain sodium hypochlorite (bleach).**
- **Do not use cleaners that contain phenol.**
- **Do not spray the selected products directly on the surfaces.**
- **All products must be used as directed by the manufacturer.**
- **Do not mix the STER 1 PLUS disinfectant with other products.**

Cleaning and disinfecting instructions.

Clean and disinfect with single-use non-abrasive paper (avoid using recycled paper) or sterile gauze. Do not use sponges or in any case, any material that can be reused.

ATTENTION!

All material used to clean and disinfect must be thrown away.

Cleaning the wheels.


Remove dust and dirt that build up inside the wheel and on the slide surface on a regular basis.

5.2. Replacing the castor wheels

Proceed as follows if one or more castor wheels on the C7-C8-C9 dentist's stool need to be replaced.

- Turn the stool upside down and then place it down with the seat resting on a table or workbench and the backrest on the outside edge of the table.
- Use a suitably-sized screwdriver to remove the castor wheel by prizing it out.
- To re-install the new castor wheel, simply insert it in the hole provided and secure it in place by banging it down hard.

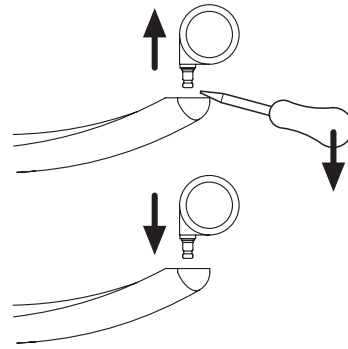5.3. Replacing the seat upholstery

Dentist's stool , model C7.

Seat: remove the seat by undoing the 4 M6x12 screws using a 4mm Allen wrench.

Backrest: remove the rear casing by prising in the slot between the casing and the upholstery (the casing is fastened with 4 snap-clips). After removing the casing, undo the 4 M6x12 screws using a 4mm Allen wrench.

Dentist's stool , model C8.

Seat: remove the seat by undoing the 4 M6x12 screws using a 4mm Allen wrench.

Dentist's stool , model C9.

Seat: remove the 2 lower casings by undoing the 4+3 M4x12 screws using a 3mm Allen wrench.

Then remove the seat by undoing the 4+4 M6x12 screws using a 4mm Allen wrench.

Backrest: remove the rear casing by prising in the slot between the casing and the vertical covering (the casing is fastened with 4 snap-clips). After removing the casing, undo the 4 M6x12 screws using a 4mm Allen wrench.

6. Trouble-shooting

P = Problem.

C = Check procedure.

A = Corrective action.

P: **The dentist's stool does not move easily on the floor.**

C: Check whether the dentist's stool wheels can all move freely or are obstructed.

A: Turn the stool upside down and clean the wheels removing any dirt that may have built up around the pivots.

NOTE: If the stool still does not move after the wheels have been cleaned, the latter need to be replaced (see paragraph 5.2).

P: **The seat cannot be lowered.**

C: Check whether the adjustment lever can move at least 3 mm up or down.

A: The handle is probably broken; call Technical Service.

P: **The seat moves down by itself without the handle being used.**

C: Check whether the adjustment lever can move at least 3 mm up or down

A: The pump is probably discharged or maintained activated by the handle; call Technical Service.

P: **The tilted backrest does not return to its initial position and/or does not provide sufficient lumbar support.**

C: The backrest turns feebly without contrast.

A: The contrast spring of the backrest is probably broken; call Technical Service.

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